

In the event you need to align the filter position with the camera, please read the following.

In some computers the mscomm32.ocx file should be registered first (run install.bat).


The Buttons 1 - 5 is the Filter Wheel control buttons.

The five up down buttons in the right is to adjust each position of the filter hole.

The calibration procedure is :

First find out which position needs to be adjusted . Then adjust the up down buttons, and click the Custom button. The Filter Wheel will flash if the data is accepted. Then click the control button to see if the position is in the right position.

The "factory" button is to set the Filter Wheel back to the factory settings.

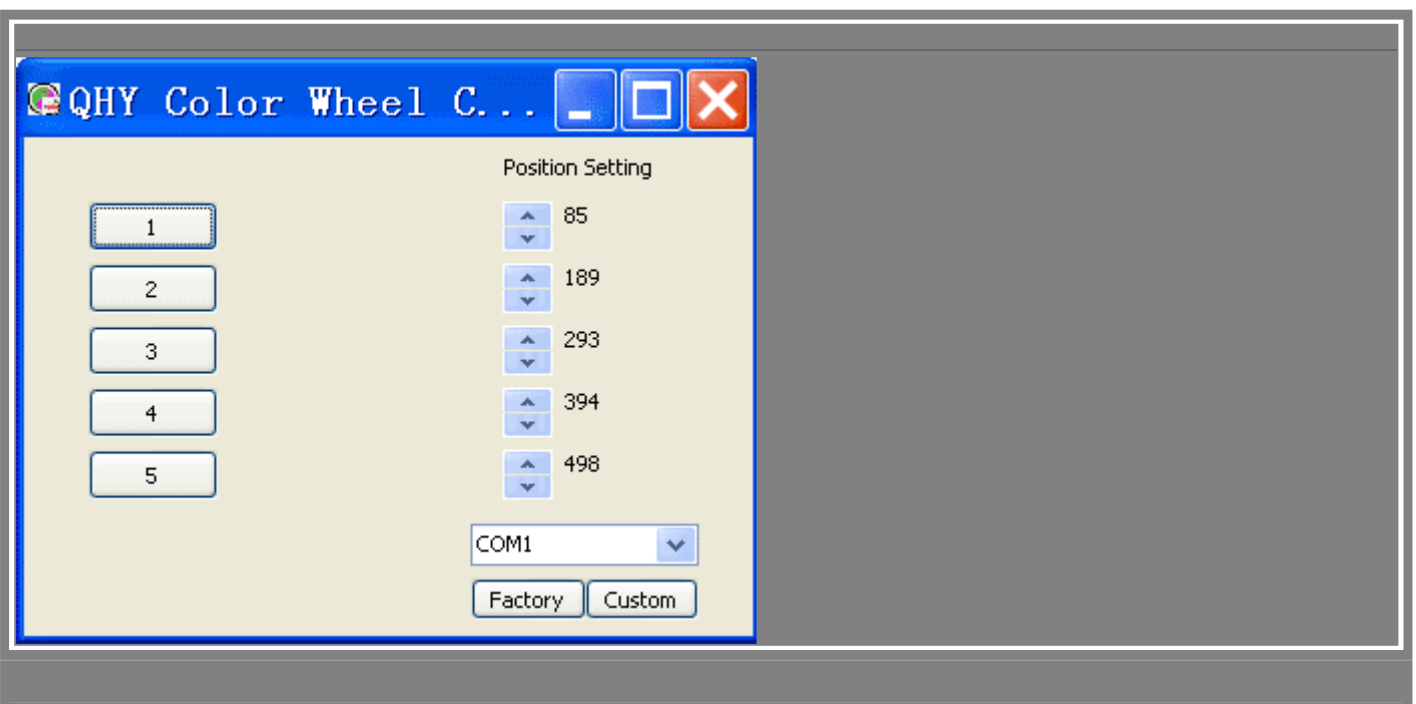

This is the control software with RS232 control.

